




Barkey *S-line*

CONVENIENT IN-LINE BLOOD & IV FLUID WARMER

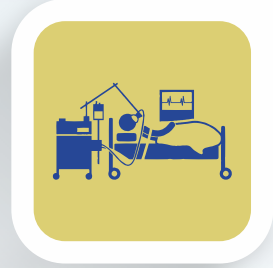

Barkey has earned its excellent reputation with exclusive engineering know-how in medical technology over the last 30 years. The proverbial quality and practical suitability of the devices and systems "made in Germany" has established itself in a wide range of applications. Real innovations and problem solutions which provide the user with specific benefits can only be developed by people who keep their ears and eyes open, who talk to specialists in practice and recognise their requirements in good time.

Barkey Sline is an innovative blood & inline IV fluid warmer, helps in preventing hypothermia in vulnerable group patients like:

- **Age:** Geriatrics (>60 Y) Infants (<6m)
- **Medical Conditions:** Hypothyroidism, strokes, parkinson's, disease, spinal cord injuries, diabetic neuropathy, long surgeries, poor nutrition, severe arthritis, trauma, severe burns, prolonged ICU patient
- **Medication:** anti-depressants, anti psychotics, sedative, drug abuse, prolonged alcoholism

Preventive temperature management has significant benefits for vulnerable groups, this reduces post-operative shivering, and recovery time



Salient Features:

- **Patented Auto-Line Technology**
Desired temperature output
- **No recurring cost of disposables**
Silicon warming sleeve to encapsulate IV drip line
- **Precise temperature management**
For safe blood transfusion
- Fluid or blood warming till the point of transfusion




Specifications subject to change without notice




Barkey S-line

CONVENIENT IN-LINE BLOOD & IV FLUID WARMER

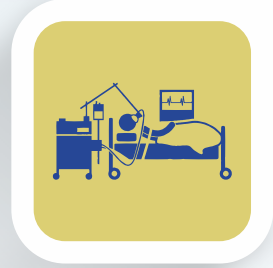

Made-to-order technology for medicine

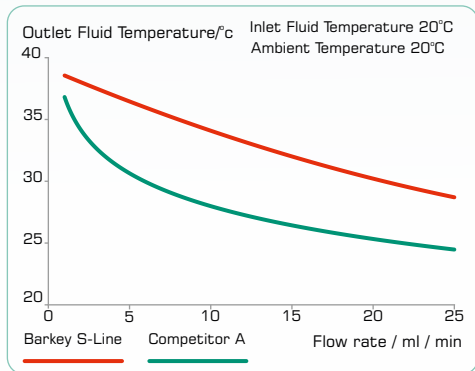

- Summary of further benefits
- Integrated temperature sensors
- Guarantee best warming and highest safety.
- Complete casing of IV-sets and therefore optimum warming
- No additional consumer articles
- Warming of blood return feeds or of dialysis fluids in combination with devices for replacement kidney therapy
- Warming of enteral and / or parenteral nutrition
- Operational safety thanks to preset target temperature value
- Display of actual temperature of the heating sleeve
- Available in a range of diameters
- Smooth surface, therefore simple to wash and disinfect
- Use of exclusively bio-compatible materials
- Optical and acoustic signal in fault situation
Exclusive innovative technological solutions at the highest quality standards and specialised benefits in practice this is the basis for the success of the holistic Barkey heating concept.

Parameter	Barkey S-line
Temperature setting	39°C set point temperature ex-works
Heating profile length	L: 1.500 mm
Internal diameters available	4.30 mm / 6.55 mm
Weight	approx. 1.300 g
Mains supply	100-240V, 50/60Hz
Connecting cable	3 m
Safety class	I protection ground Degree of protection BF external and internal application to patients with insulated application parts
Humidity protection	IP X2
Certification	CE 0123
Country of origin	Germany



Specifications subject to change without notice